# Michelson Experiment Summary

Yongao Hu

MIT Department of Physics

March 3, 2025

#### Outline

- The History of Michelson Interferometer
- Theoretical background: Interference of Light
- Experimental Setup: the Interferometer
- Oata Analysis: Detection of Minima and Maxima
- Sesults: Visibility and Wavelength of the Laser
- Conclusion

#### The History of Michelson Interferometer

- The Michelson interferometer is designed by Albert Michelson.
- Historically, it disproved the existence of the "luminiferous aether" and contributed to the discovery of relativity.
- Interferometers are used in:
  - Precision measurements
  - LIGO for detecting gravitational waves [1]

#### Interference of Light Waves

• The interference pattern depends on the phase difference:

$$\Delta \phi = 2\pi \frac{\Delta x}{\lambda}$$

• Intensity of interference:

$$I = I_0 \cos^2(\Delta \phi)$$

- Maximum interference:  $\Delta \phi = 2\pi n$
- Minimum interference:  $\Delta \phi = (2n+1)\pi$

## Interference of Light Waves

 The interference pattern depends on the phase difference:


$$\Delta \phi = 2\pi \frac{\Delta x}{\lambda}$$

• Intensity of interference:

$$I = I_0 \cos^2(\Delta \phi)$$

- Maximum interference:  $\Delta \phi = 2\pi n$
- Minimum interference:  $\Delta \phi = (2n+1)\pi$

#### Constructive Interference

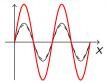

#### Destructive Interference

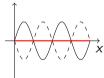


Figure: Solid/dash lines: individual waves. Red lines: resultant waves

#### Visibility of the Interference Pattern Measures how Resolved the Pattern is

• Visibility measures how well the pattern is resolved:

$$V = \frac{I_{\text{max}} - I_{\text{min}}}{I_{\text{max}} + I_{\text{min}}}$$

- Ideal case: V=1
- Low visibility ⇒ poor alignment or coherence issues

#### Piezoelectric Transducer Controls the Path Difference

- A piezoelectric transducer moves a mirror in response to voltage.
- This changes the path length  $\Delta x$  proportional to voltage.
- The wavelength of light is found using:

$$\lambda = 4\Delta V \times \frac{\Delta L}{V}$$

## Michelson Interferometer Setup

- Light source: laser
- 50/50 beam splitter splits the light into two paths.
- One mirror is stationary, the other mounted on a piezoelectric transducer.
- Detector: photodiode (voltage output proportional to intensity).
- Signal generator applies a triangular voltage to move the mirror.

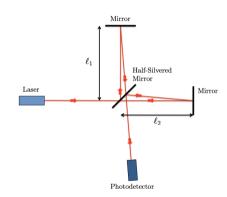


Figure: Schematic of the Michelson interferometer. Adapted from [2].



## Data Processing and Interference Pattern Visibility

- Voltage signals recorded from signal generator and photodiode.
- Used findpeak in SciPy to locate primary maxima and minima.
- Applied KMeans clustering (SciKit-Learn) to group voltage values.

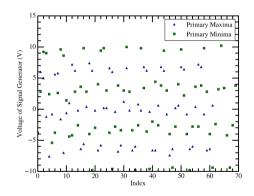


Figure: Voltage at primary maxima and minima.

# Visibility $0.179 \pm 0.001$ Not Perfect but Significant

- ullet Maximum voltage:  $(1.592 \pm 0.002)$  V
- ullet Minimum voltage:  $(1.110 \pm 0.002)$  V
- Visibility:

$$V = 0.179 \pm 0.001$$

Imperfect alignment likely lowered visibility.

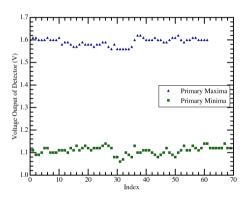


Figure: Voltage output of photodiode corresponding to maxima and minima



## Wavelength $601 \pm 37 (\text{sys}) \pm 10 (\text{stat}) \text{nm}$ Matches the Color of the Laser

- Measured voltage difference between peaks:  $(3.192 \pm 0.055) \text{ V}$
- Calculated wavelength:

$$\lambda = 601 \pm 37 ( ext{sys}) \pm 10 ( ext{stat}) \; ext{nm}$$

 Matches expected wavelength of the orange laser.

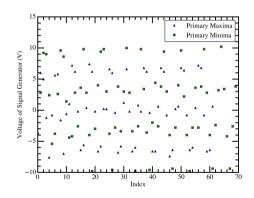


Figure: Voltage at primary maxima and minima.

#### Uncertainty Dominated by Calibration of the Piezoelectric

• Systematic error: Calibration of the piezoelectric device

$$\delta\lambda=$$
 37 nm

• Statistical error: Peak detection noise

$$\delta\lambda=10\;\mathrm{nm}$$

#### Conclusion and Future Improvements

• Successfully measured the wavelength of an orange laser:

$$601 \pm 37 ( ext{sys}) \pm 10 ( ext{stat}) \ ext{nm}$$

- Demonstrated interference and wave nature of light.
- Future improvements:
  - Improve mirror alignment for higher visibility.
  - Use a more precise piezoelectric transducer.
  - Implement automated feedback to stabilize fringes.



#### Acknowledgments

Special thanks to J. Lewis, M. Bohdan, and V. Tran for assistance. Thanks to the MIT 8.13 teaching team for guidance.

#### References I



B. P. Abbott, et al. (LIGO Scientific Collaboration, and Virgo Collaboration). Observation of gravitational waves from a binary black hole merger. *Physical Review Letters*, 116:061102, 2016.



MIT Department of Physics.

Optical interferometry.

Technical report, MIT Department of Physics, August 2023.

